
BBS300

Personenwaage

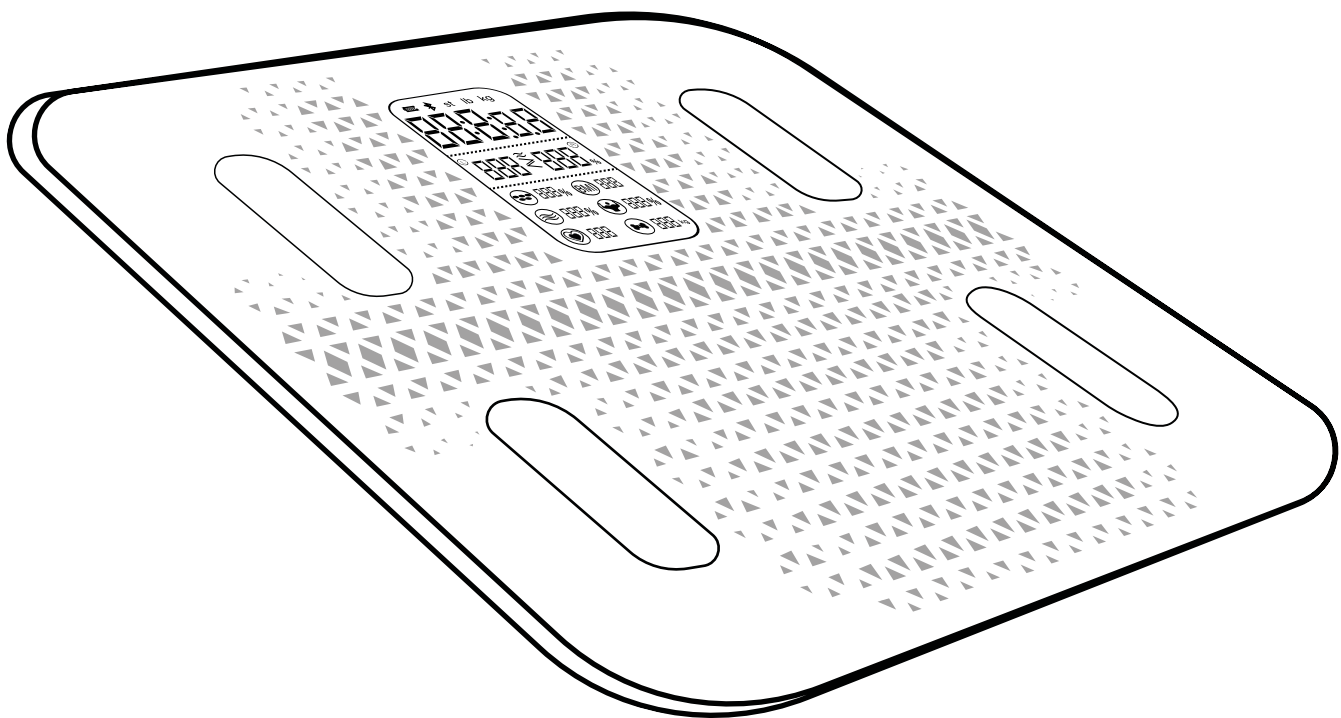DE

EN

ES

FR

IT

NL

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DEUTSCH

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ENGLISH

Dear customer,

We are pleased you chose a device from the SPORTSTECH product range. SPORTSTECH sports equipment offers you high quality and new technology.

In order to fully use the potential of your device and be able to enjoy it for many years, please read this manual carefully before starting up and beginning of training, and use the device according to the instructions. The operational safety and function of the device can only be guaranteed if the safety instructions in this user manual are observed. We shall not assume any liability for damages resulting from improper use or incorrect operation.

Please ensure that all persons using the device have read and understood the user manual.

Keep the user manual in a safe place to be able to access it at any time if needed.

Our **video tutorials** for you!

Assembly, use, disassembly.

1. Scan the QR code
2. Watch videos
3. Start fast and safely


Link to the videos:

https://service.innovamaxx.de/bbs300_video

We are also on social media!

Get the latest product information, training content and much more on our:

Instagram-page

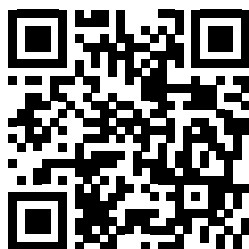Facebook-page

<https://www.instagram.com/sportstech.de>

<https://www.facebook.com/sportstech.de>

INDEX

1. INTRODUCTION.....	24
2. IMPORTANT SAFETY INSTRUCTIONS.....	26
3. OVERVIEW	28
4. DISPLAY AND FUNCTIONS	29
5. INSTALLATION	30
6. BODY ANALYSIS	31
7. NORMAL WEIGHING FUNCTION	32
8. APPLICATION SOFTWARE (APP).....	33
9. RESULT DISPLAYS	34
10. TROUBLESHOOTING.....	35
11. TECHNICAL DATA.....	35
12. DISPOSAL	36
13. CONFORMITY	37

1. INTRODUCTION

The SPORTSTECH body analysis scale BBS300 measures using the Bioelectric Impedance Analysis (short BIA) muscle mass, bone mass, body fat and body water content and gives you the percentage distribution of these values in relation to your total body weight.

Based on the evaluation of this personal data, the scale also shows you your body mass index (BMI) and your heart rate. The device thus offers you a reliable and precise overview of your body composition.

Intended Use

This product is intended for private use and is not suitable for medical, therapeutic or commercial use.

The maximum load capacity of 180 kg must not be exceeded.

Attention! This device is not suitable for people with medical implants such as a pacemaker or other electronic medical devices. Their function could be impaired.

How does the BIA measurement work with the SPORTSTECH BBS300?

The advanced BIA measurement takes place via a weak and harmless electrical pulse, which is sent through the body by the 4 integrated electrodes within just a few seconds. The electrical resistance in the body is measured, whereby the fat tissue and the bones have more resistance and therefore less electrical conductivity than the muscle tissue and body water.

The person-dependent, individual factors such as gender, height and age are taken into account in order to determine your personal body composition as accurately as possible.

Please note that the determined values of the SPORTSTECH bathroom scale only represent approximate values to the real, medical analysis values of the body. The exact determination of muscle percentage, bone mass, body fat and body water can only be carried out using the methodology of a medical specialist.

13 different measurements and values

in the "Fitdays" app

- **Weight**
- **BMI (Body Mass Index)**
- **Body fat percentage**
- **Body weight without fat**
- **Subcutaneous fat**
- **Visceral fat**
- **Body water**
- **Muscle percentage**
- **Muscle mass**
- **Bone mass**
- **Protein balance**
- **Basal metabolic rate (BMR)**
- **Metabolic age**

EN

2. IMPORTANT SAFETY INSTRUCTIONS

These basic precautions should always be followed, including the following safety instructions when using this device. Read all instructions before using this equipment.

1. It is important to read this entire manual before using the equipment. Safe and effective use can only be achieved if the equipment is maintained and used properly. It is your responsibility to ensure that all users of the equipment are informed of all warnings and precautions.
2. Keep children away from the packaging material of the product, there is a risk of suffocation.
3. Attention! This device is not suitable for people with medical implants such as a pacemaker or other electronic medical devices. Their function could be impaired.
4. Do not use the device during pregnancy or if you think you may be pregnant.
5. Use the device on a firm, level and non-slippery surface.
6. Only step on the intended stand area and not on the actual display.
7. Caution! Danger of Slipping! Do not use the device if its surface is damp or your feet are wet.
8. Caution! Risk of Falling! Do not stand on one side on the outer edge of the scale.
9. Make sure that small children and pets are kept away from the device to avoid injury.
10. The individual components of this device must not be disassembled, converted or repaired.
11. Please note that there may be measurement variables due to technical reasons (0-60kg: $\pm 300g$, 60-120kg: $\pm 400g$, 120-180kg: $\pm 500g$), as the device is not a scale that has been calibrated for professional medical use.
12. The safety of the equipment can only be maintained if it is regularly examined for damage and/or wear and tear. Repairs may only be carried out by specialist personnel.
13. Always use the equipment as indicated. If you find any defective components while assembling or checking the equipment, or if you hear any unusual noises coming from the equipment during use, stop immediately. Do not use the equipment until the problem has been rectified.
14. This product is intended for private use and is not suitable for medical, therapeutic or commercial use.

NOTES ON BATTERY SAFETY

1. Protect the bathroom scale from excessive sunlight, heat, fire, etc.: Otherwise the built-in battery could explode.
2. Please never try to open the device yourself. If the bathroom scale may be open due to damage, do not touch the contents, especially the battery poles!
3. Please do not open the body of the bathroom scale to change the battery yourself. The battery may only be changed by specialist personnel from Sportstech Brands Holding.
4. To recharge the battery, please use only the charger supplied. If you perceive a smell, or if the battery heats up, please stop the charger immediately.
5. The batteries cannot be removed or replaced. They are fixed to the device.
6. Do not charge the device overnight or unattended.

NOTES ON CHARGER SAFETY

7. The battery charger must not be used by children.
8. Children are not allowed to play with the charger.
9. Children must not clean or service the device. If the charger cable is defective, it must be replaced by the manufacturer, a specialist dealer or another qualified person. Otherwise there is danger.
10. It is crucial to protect electrical parts from moisture. When cleaning the appliance, or when it is in use, make sure that no liquid is dripping or flowing into the appliance. Otherwise, you may suffer an electric shock. Do not keep the charger under water. Read the instructions in the section on cleaning and care thoroughly..
11. When using the charger, be sure to periodically check for damage. You can only use a defective charger again if it has been refurbished.
12. Do not use the charger if it has taken a hit or a knock: for example, when it fell or is damaged for another reason.
13. If the bathroom scale is defective, do not connect a charger. There is a risk of injury from electric shock.
14. Do not disassemble the charger under any circumstances..
15. Do not use the charger near explosive or flammable objects. Otherwise, there is risk of explosions and fire.
16. Please charge the battery only in a confined space, as the charger is designed for indoor use only. Otherwise, an electric shock might occur.
17. Please remember that the charger and the charger port must be connected correctly.
18. Always make sure that the charger port is in fact free. Wetness and dirt may damage the charger port.
19. Please also read our section „Disposal“.

3. OVERVIEW

SCOPE OF DELIVERY:

- Personal scale BBS300 with Li-polymer battery cell
- Micro USB cable (60 cm long)
- User Guide

General information on the weighing process with the SPORTSTECH BBS300

- For the most accurate measurement possible, stand still during the measurement process and do not move if possible. Do not support yourself or hold on to anything. Your feet should be positioned centrally and symmetrically on the electrodes of the scale.
- If possible, always carry out the weighing process at the same time of day and under the same conditions (undressed, sober) in order to obtain comparable results.
- The body analysis values can only be determined barefoot. Since completely dry soles of the feet have too low a conductivity, you can moisten them very weakly if necessary.
- The measured values in a body analysis with bioimpedance measurement methods depend on several factors and can vary considerably over the course of the day. Weighing in the evening (between 6 and 8 p.m.) is, therefore, recommended as an optimum time for a body analysis measurement. In the morning, the body is drained, which means that the body resistance and thus the measured body fat percentage is higher, but the muscle percentage is accordingly lower. During the day, the body's water balance returns to a normal level. You can find more detailed information in section 9 Result Displays.


Notes on determined measurement results

The results can be falsified if the device

- is on an uneven surface.
- is used in different locations.
- is close to mobile phones.
- is placed in locations with high humidity.
- is close to heating sources or in the range of fluctuating temperatures.
- is exposed to shocks or vibrations.

Avoid direct sunlight.

4. DISPLAY AND FUNCTIONS

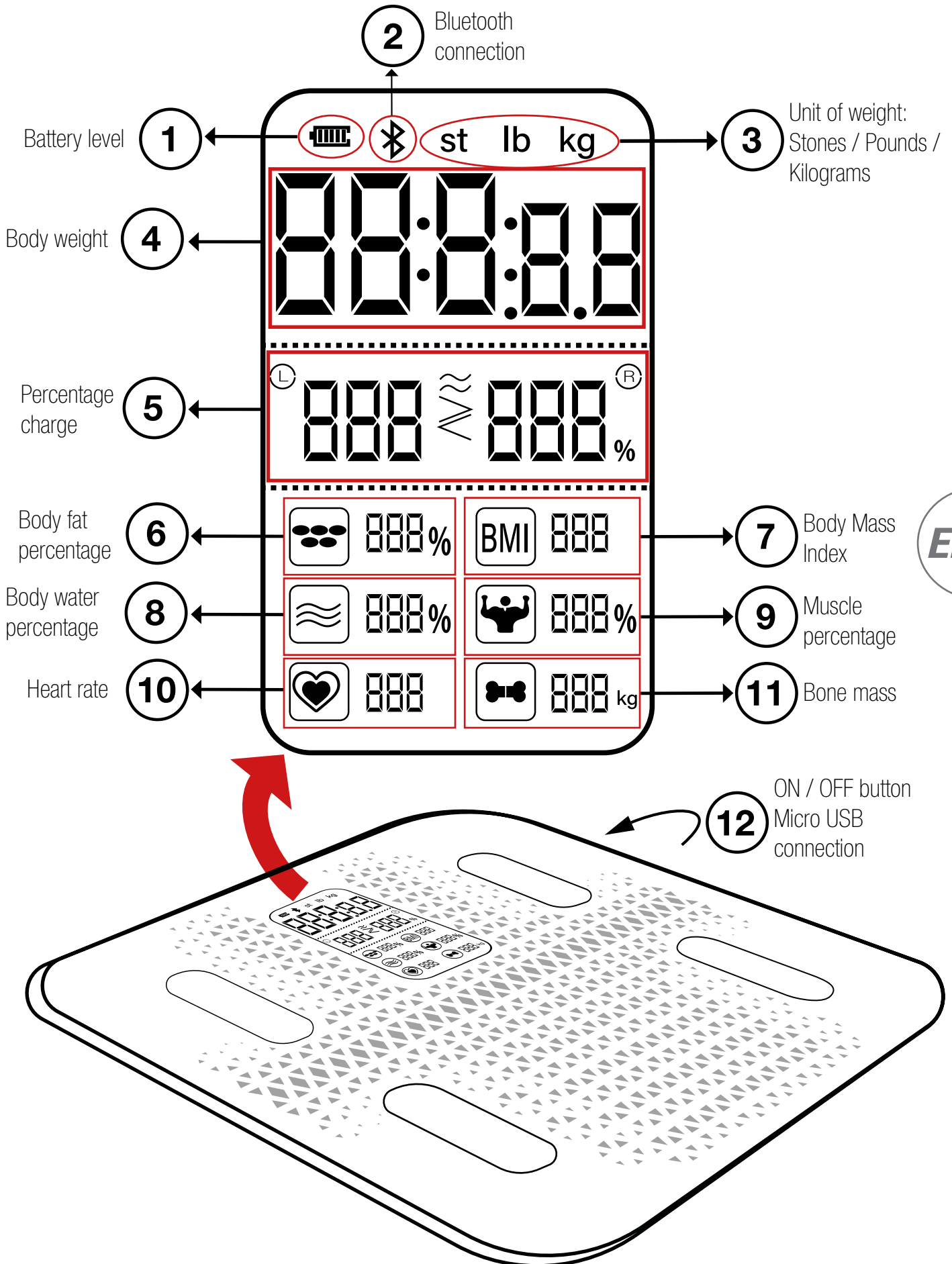5. INSTALLATION

1. Charging the battery

Before using for the first time, fully charge the integrated battery using the micro USB cable.

The battery symbol on the top left of the display shows you the current charge level. With average use, the battery life lasts at least 30 days.

2. Personal Settings


The body analysis measurement can display the determined values very precisely and individually based on your personal settings in the "Fitdays" app. The app offers you 13 different measurement values to keep an eye on during your fitness and health regimes and to display any progress and changes for you.

Before you carry out the first measurement with your SPORTSTECH balance, please install the free "Fitdays" app on your smartphone and create a user profile. After the scale has been connected to the app once via Bluetooth, the scale shows you the analysis values on its display. The "Fitdays" app offers you further diverse measurement values, individual setting options and long-term comparison values from your measurements. You can also select your preferred weight unit via the app.

You can find more information on this under point 8 Application Software (App).

3. Cleaning and Care

Please only use a dry, soft cloth to clean the scale. You can also use a cloth that has been lightly moistened with water or a mild detergent. Do not hold the device under running water. Refrain from using aggressive cleaning agents such as alcohol, abrasives or solvents.

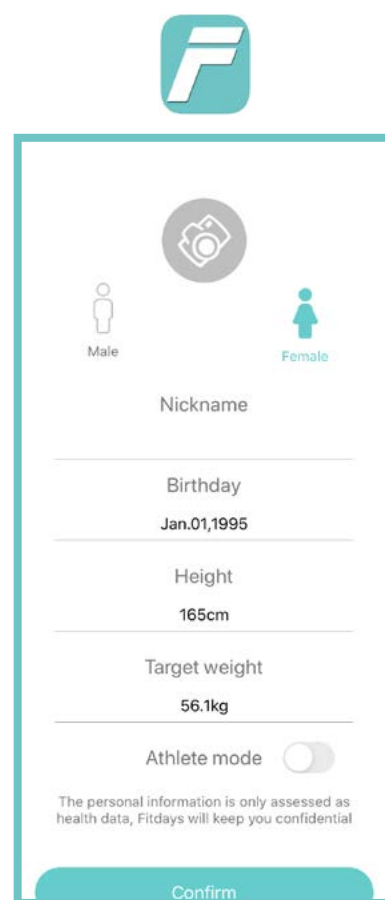6. BODY ANALYSIS

Start a measurement

1. Open the "Fitdays" app on your smartphone and make sure that the SPORTSTECH scale is connected via Bluetooth.
2. Make sure that the scale is on a level, firm surface. The body analysis function is only available when you stand barefoot on the scale.
3. Press the ON / OFF button once on the right side of the device (next to the micro USB connection) so that the display lights up.
4. Mount the scales barefoot and position your feet symmetrically on the 4 silver electrodes. The body analysis starts automatically.

Factors that can influence a body analysis

Basically, the measured values of each person fluctuate during the day according to the distribution of body water. In women, hormonal fluctuations can change the water balance in the body, which can lead to deviations in the measurement results. It is not recommended to measure yourself immediately after strenuous activities, sauna sessions, bathing, alcohol consumption and eating and drinking, as this can influence the values to be determined.

The following persons or groups of people could receive different or compromised results when determining the body analysis values:

- Competitive athletes and bodybuilders
- Pregnant women
- Children under 10 years
- People on dialysis treatment, with fever, edema symptoms or osteoporosis
- People who are taking vasoconstrictor or vasodilator medication
- People taking cardiovascular medication
- People whose overall body size is influenced by considerable anatomical deviations in the legs (difference in leg length)

7. NORMAL WEIGHING FUNCTION

Normal weighing function (weight measurement only)

It is also possible to use the device for normal weight measurement. To do this, press the ON / OFF button once on the right side of the device (next to the micro USB port) so that the display lights up.

Now stand on the scales with stockings or shoes and you will only be shown your body weight.

8. APPLICATION SOFTWARE (APP)

You can download the "Fitdays" app using the QR codes shown here or search in your store. Minimum compatibility specifications for mobile device; Android 5.0 or higher with Bluetooth 4.0, iOS 8.0 or higher.

Android

Fitdays

iOS

Notes on using the App:

You should only install apps from sources you trust. Sportstech Brands Holding cannot accept any liability for software provided by third parties.

Please remember to turn on the Bluetooth feature on your device.

Android & iOS

1. After downloading the app, register in just a few steps.

2. The app now prompts you to establish a Bluetooth connection and guides you through the weighing process after successful connection is established.

9. RESULT DISPLAYS

Body fat percentage

Age (in years)		Women			Men	
	WELL	MEDIUM	BAD	WELL	MEDIUM	BAD
< 20	17-22 %	22-27 %	> 27 %	12-17 %	17-22 %	> 22 %
20-30	18-23 %	23-28 %	> 28 %	13-18 %	18-23 %	> 23 %
30-40	19-24 %	24-29 %	> 29 %	14-19 %	19-24 %	> 24 %
40-50	20-25 %	25-30 %	> 30 %	15-20 %	20-25 %	> 25 %
> 50	21-26 %	26-31 %	> 31 %	16-21 %	21-26 %	> 26 %

Body water percentage

Age (in years)		Women			Men	
	VERY WELL	WELL	BAD	VERY WELL	WELL	BAD
10-100	> 60 %	45-60 %	< 45 %	> 65 %	50-65 %	< 50 %

Muscle percentage

Age (in years)		Women			Men	
	LOW	NORMAL	HIGH	LOW	NORMAL	HIGH
15-19	< 35 %	35-41 %	> 41 %	< 43 %	43-56 %	> 56 %
20-29	< 34 %	34-39 %	> 39 %	< 42 %	42-54 %	> 54 %
30-39	< 33 %	33-38 %	> 38 %	< 41 %	41-52 %	> 52 %
40-49	< 31 %	31-36 %	> 36 %	< 40 %	40-50 %	> 50 %
50-59	< 29 %	29-34 %	> 34 %	< 39 %	39-48 %	> 48 %
60-69	< 28 %	28-33 %	> 33 %	< 38 %	38-47 %	> 47 %
> 70	< 27 %	27-32 %	> 32 %	< 37 %	37-46 %	> 46 %

Bone mass in kilograms

		Women			Men	
Weight (in kg)	< 45	45-60	> 60	< 60	60-75	> 75
Bone mass (in kg)	1,8	2,2	2,5	2,5	2,9	3,2

* The values shown are standardized and therefore not generally applicable. Deviations and anomalies are possible in your individual measurements, so always contact a doctor for medical advice.

10. TROUBLESHOOTING

Problem	Possible Cause	Solution
Measurement results are abnormal: <ul style="list-style-type: none"> • Too low • Too high • Large difference between 2 real-time measurements 	<ul style="list-style-type: none"> • The device is on an uneven surface • Feet are too dry (poor conductivity) • Wrong posture or support 	<ul style="list-style-type: none"> • Place the scale on a firm, level surface. • Moisten your feet mildly with a damp cloth. • Step on the scales barefoot and ensure the same foot position on both sides.
The LED display does not turn on after you press the ON / OFF button.	The battery is empty.	Charge the battery using the micro USB cable.
The scales do not perform a body analysis.	Weighing is done wearing socks or shoes.	Mount the scale barefoot and keep your feet in contact with the 4 electrodes during the measurement.

11. TECHNICAL DATA

Product model	BBS300
Dimensions	300 x 300 x 25 mm
Display dimensions	62 x 109 mm
Units of weight	Pound / Stone / Kilogram
Measuring range	5kg-180kg / 11lb-397lb / Ost 11lb-28st 5 lb
Net weight	2,16 kg
Maximum load capacity	180 kg
Battery	Model: JY502030 Type: Li-polymer rechargeable battery cell Power: 300 mAh Voltage: 3,7 V Manufacturer: Zhongshan Jingyu Electronic Technology Co.,Ltd.
Battery charger	Micro USB cable
Radio frequency	2400-2483.5MHz
Maximum transmission power	≤20 dbm

12. DISPOSAL

Dear customers,

As an end user, you are legally obliged to collect waste electrical and electronic equipment and batteries as well as accumulators separately from municipal waste. We want to make this as easy as possible for you. Please note the information on this page and on the back.

Notes on EU directive 2012/19/EU (in Germany implemented in ElektroG (Law about electric and electronic devices); elsewhere implementation may differ)

We advise owners of electric and electronic equipment that waste electrical and electronic equipment must be collected separately from municipal waste in accordance with the applicable legal regulations. The symbol of a crossed-out garbage can, which is shown on waste electrical and electronic equipment, also indicates the obligation to collect this waste separately.

B. Notes on disposal and symbols for an electrical appliance based on EU Directive 2006/66/EU (in Germany in accordance with the Battery Act (BattG); elsewhere, implementation may differ)

Old batteries and accumulators must not be disposed of with household waste. As an end user, you are legally obliged to return used batteries. The symbol with the crossed-out garbage can means that you are not allowed to dispose of old batteries in the household waste.

Use the official collection points for this, such as the collection points of public waste disposal carriers. You can also return used batteries to retail outlets, provided they sell batteries. If the symbol Cd, Hg, or Pb is displayed under the symbol of the crossed-out garbage can, this is an indication that the battery contains either cadmium, mercury or lead. These are toxic heavy metals that not only pose significant risks to the environment but also to human health.

Recycling loop

Packaging material can be conveyed back to the raw material cycle. Dispose of packaging material in accordance with legal provisions. Information can be retrieved from the return or collections systems of your community.

Information about included batteries


Information according to § 4 paragraph 4 Elektrogesetz

Battery type:	Li-ion battery
Chemical system:	Li-NMC lithium nickel cobalt manganese oxide
Battery manufacturer	Dongguan Lean Power New Energy Technology Co.,Ltd

Information on the safe removal of the batteries or accumulators

Tools needed: Phillips screwdriver, flathead screwdriver

- The battery/rechargeable battery is located under a cover of the device and is also connected to the main circuit board with a cable. The removal is not possible without a certain effort.
- If you can remove the battery/rechargeable battery professionally, dispose of the battery/rechargeable battery and, if applicable, the device separately at the appropriate official collection point.
- If you cannot remove the battery/rechargeable battery, hand in the entire device to an official collection point. Please also note the locally applicable regulations.

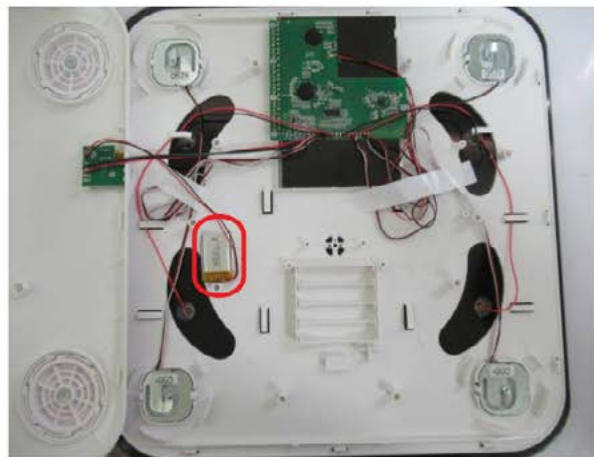13. CONFORMITY

Sportstech Brands Holding hereby declares that the radio system type is compliant with EU Directive 2014/53 / EU. The full text of the EU declaration of conformity is available at the following internet address: https://service.innovamaxx.de/bbs300_conformity

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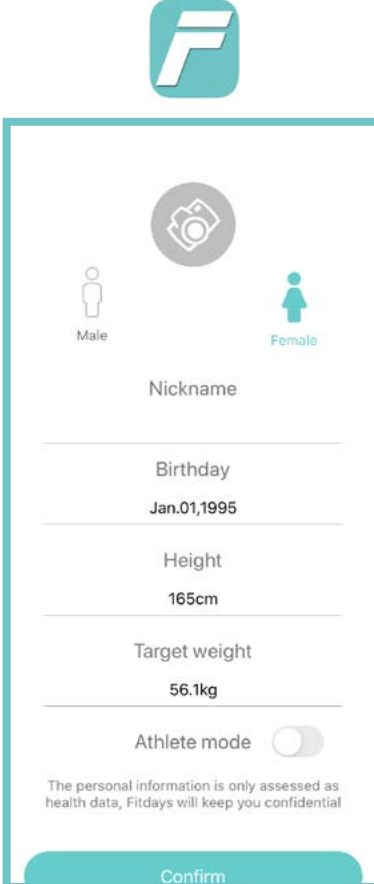
The image shows a screenshot of the Fitdays app's user profile setup screen. At the top is the Fitdays logo, a stylized 'F' in a teal square. Below it is a camera icon for profile picture. There are two gender selection options: 'Male' with a grey icon and 'Female' with a teal icon. The 'Female' option is selected. Below these are fields for 'Nickname', 'Birthday' (set to Jan.01,1995), 'Height' (set to 165cm), and 'Target weight' (set to 56.1kg). There is an 'Athlete mode' toggle switch which is currently turned off. At the bottom, there is a disclaimer: 'The personal information is only assessed as health data, Fitdays will keep you confidential' and a teal 'Confirm' button.

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Sportstech Brands Holding GmbH