
Relaxknow

# RXM50

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Massager

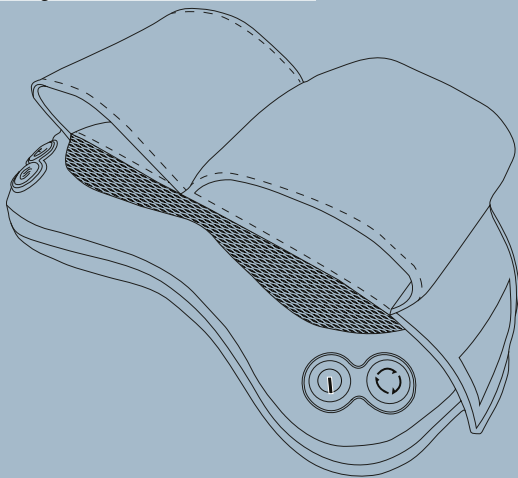DE

EN

ES

FR

IT

NL

# DEUTSCH

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# ENGLISH

Safety Instructions.....	12
Technical Specifications and Delivery .....	13
Warning .....	14
Function.....	14
Operation .....	15
Trouble-shooting.....	15
Warranty .....	16
Cleaning and Care .....	17
Disposal .....	18
Declaration of Conformity.....	51

This manual represents a status at press date. It is subject to changes without notice and to errors and omissions.

## Our video tutorials for you! Assembly, use, disassembly.

1. Scan the QR code
2. Watch videos
3. Start fast and safely

Link to the videos:

[https://service.innovamaxx.de/rxm50\\_video](https://service.innovamaxx.de/rxm50_video)

## We are also on social media!

Get the latest product information, training content and much more on our:

Instagram-page

<https://www.instagram.com/sportstech.de>

Facebook-page

<https://www.facebook.com/sportstech.de>

## SAFETY INSTRUCTIONS

Please read these operating instructions carefully before using the device. Use the RXM50 Massager only as described in this manual. Any other use is deemed incorrect and may result in damage to property or personal injury. The RXM50 Massager is not a children's toy.

The RXM50 Massager is designed exclusively for massaging body areas. Store your RXM50 Massager in a clean and dry place protected from major temperature changes and other potentially harmful factors. Additionally, the RXM50 Massager is designed for indoor use only.

- Do not use the RXM50 Massager when it is wet or damp.
- Dismantling or rebuilding the RXM50 Massager may cause damage and injury
- Use the RXM50 Massager for approx. 15 minutes and then switch it off for approx. 20 minutes.
- If the power adapter or the adapter is damaged, the RXM50 Massager must not be used!

*The manufacturer or dealer assumes no liability for damage caused by improper or incorrect use.*

**Our brands** - Discover many other interesting products:

**Sportstech  
Fitness**

<https://www.sports-tech.uk.com/>

**Bluewheel  
Electromobility**

<https://www.bluewheel.de/>

**Deskfit**

<https://www.deskfit.de/>

**Relaxxnow**

via Sportstech

## TECHNICAL SPECIFICATIONS AND DELIVERY

Dimensions: 310 x 190 x 90 mm

Net weight: 1,0 Kg

Massager RXM50

AC adapter

Vehicle power adapter

Operation manual

### Technical data of the AC adapter

Manufacturer's name or trade mark, trade register number and address	Shenzhen Zhongli Power Technology Co. Ltd Building F, Lihao Yongda High-tech Park, Gongming Street, Guangming New District, Shenzhen, China	
Model Identifier	ZL-024WL1202000EU01	
	<b>Nominal or measured value</b>	<b>Unit</b>
Input Voltage	100..240	V AC
Input AC Frequency	50..60	Hz
Output Voltage	12	V DC
Output Current	2	A
Output Power	24	W
Average active efficiency	89.03	%
Efficiency at light load (10%)	86.21	%
No-load power consumption	0.039	W

## **WARNING**

Do not insert any body parts into the device to avoid injuries and damage!

The following people should not use the RXM50 Massager:

- People with osteoporosis
- People with a pacemaker
- People with heart disease
- Pregnant women and people with wounds and/or injuries
- And children under 12 years


## **FUNCTION**

The RXM50 Massager is deemed suitable for the following people:

- People sitting / standing in the same position for a long time (for example sitting at a desk in the office for a long time)
- Nervous and / or tense people
- People who suffer from sleep disorders
- People with a stiff neck or stiff shoulder parts
- People with headaches or dizziness
- People with 'sticky' fascia

### **Before plugging into the main electricity + battery mode**

1. Remove the packaging material and all protective films.
2. Unwind the cable of the power plug completely to straighten it.
3. Insert the connector into the hollow plug of the power supply.
4. Plug the power plug into a properly installed power outlet.

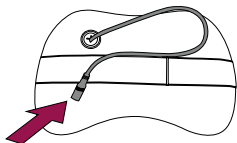You can use the RXM50 Massager as shown for the whole body:

NECK

WAIST

ABDOM


LEG

ARM

CALF

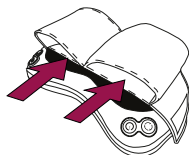Loop for hands or feet

## OPERATION

Power: turns device on/off

Direction: changes direction

Heat: turns on/off heating

Speed: slow / middle / fast

## ! TROUBLE-SHOOTING

- If the RXM50 Massager can not be switched on, please check whether the plug is correctly inserted and check the power supply.
- If the sound of the RXM50 Massager gets progressively louder, it may be because it is in continuous use for too long. Switch off the RXM50 Massager for at least 30 minutes before restarting.
- If the RXM50 Massager stops during use, it may be because the engine has overheated. Leave the RXM50 Massager turned off for at least 30 minutes before restarting.

Problem	Possible cause	Problem Resolution
The RXM50 Massager makes strange or loud noises.	The massage balls are electrically operated such that their activation is audible.	No correction necessary.
	The noise becomes increasingly louder.	Switch off the RXM50 Massager for at least 30 minutes.
	The massage balls rub on the fabric during normal operation.	No correction necessary.
The RXM50 Massager does not work.	The connectors are not connected properly.	Check the connectivity of the power connector and the power supply connector.

**IMPORTANT:** In the event of a fault, do not attempt any repairs yourself. There are no warranty claims in such cases.

## General information

The RXM50 Massager is intended for private use only and not for medical or commercial use. It is suitable for full body massage: neck, shoulder, back, legs, stomach, feet and arms. If you have any health concerns, please talk to your doctor before using the RXM50 Massager.

The massage should always be a pleasant experience. If you find the massage uncomfortable or causing pain, please change the intensity of the massage and/or adjust the tension of the arm loops to achieve a satisfactory level of comfort.

Do not use the RXM50 Massager while driving.

**After 15 minutes of uninterrupted massage, the device switches off automatically.**

## CLEANING AND CARE

Before cleaning the RXM50 Massager, be sure to fully turn off the power and unplug the adapter.

Please note:

- Note that the RXM50 Massager is not washable per se!
- Always use a dry cloth for cleaning.
- Under no circumstances should any liquid get into the device or the accessories, otherwise a short circuit may occur.
- Minor, small stains may be removed only with a slightly damp cloth or sponge (with possibly a little soapy water).
- Do not use aggressive or solvent based cleaners. Do not use brushes with metal or nylon bristles, as well as sharp or metallic cleaning items such as knives and similar sharp objects.
- After cleaning, do not use the RXM50 Massager until it has completely dried.
- Keep the RXM50 Massager in a protective cover when not in use to prevent the accumulation of dust.

## DISPOSAL

### Dear customers,

As an end user, you are legally obliged to collect waste electrical and electronic equipment and batteries as well as accumulators separately from municipal waste. We want to make this as easy as possible for you. Please note the information on this page and on the back.

**Notes on EU directive 2012/19/EU** (in Germany implemented in ElektroG (Law about electric and electronic devices); elsewhere implementation may differ)

We advise owners of electric and electronic equipment that waste electrical and electronic equipment must be collected separately from municipal waste in accordance with the applicable legal regulations. The symbol of a crossed-out garbage can, which is shown on waste electrical and electronic equipment, also indicates the obligation to collect this waste separately.

**B. Notes on disposal and symbols for an electrical appliance based on EU Directive 2006/66/EU** (in Germany in accordance with the Battery Act (BattG); elsewhere, implementation may differ)

Old batteries and accumulators must not be disposed of with household waste. As an end user, you are legally obliged to return used batteries. The symbol with the crossed-out garbage can means that you are not allowed to dispose of old batteries in the household waste.

Use the official collection points for this, such as the collection points of public waste disposal carriers. You can also return used batteries to retail outlets, provided they sell batteries. If the symbol Cd, Hg, or Pb is displayed under the symbol of the crossed-out garbage can, this is an indication that the battery contains either cadmium, mercury or lead. These are toxic heavy metals that not only pose significant risks to the environment but also to human health.

### Recycling loop

Packaging material can be conveyed back to the raw material cycle. Dispose of packaging material in accordance with legal provisions. Information can be retrieved from the return or collections systems of your community.

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Unterszeichnet für und im Namen der Sportstech Brands Holding GmbH von: Ali Ahmad
**Sporttech Brands Holding GmbH**