
Relaxknow

RXM200

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Neck and Shoulder Massager

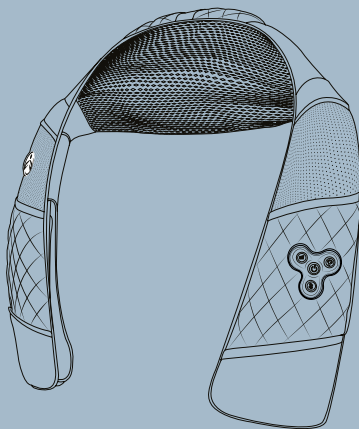DE

EN

ES

FR

IT

NL

DEUTSCH

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
ENGLISH

Technical specifications and delivery	13
Safety instructions	13
Warning	14
Function	14
Service	15
Trouble-Shooting	16
Warranty	17
Disposal	19
Declaration of conformity	20

Our **video tutorials** for you!

Assembly, use, disassembly.

1. Scan the QR code
2. Watch videos
3. Start fast and safely

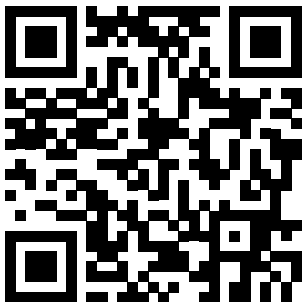Link to the videos:

https://service.innovamaxx.de/rxm200_video

We are also on social media!

Get the latest product information, training content and much more on our:

Instagram-page

Facebook-page

<https://www.instagram.com/sporttech.de>

<https://www.facebook.com/sporttech.de>

TECHNICAL SPECIFICATIONS AND DELIVERY

Power: 30 W

Voltage: AC100-240V 50 / 60Hz - DC12 V, 2.5 A

Dimensions: 400 x 450 x 200mm

Net weight: 1,9 Kg

Massager RXM200

Vehicle AC adapter

AC adapter

Operation manual

SAFETY INSTRUCTIONS

Please read these operating instructions carefully before using the device. Use the RXM200 Massager only as described in this manual. Any other use is deemed incorrect and may result in damage to property or personal injury. The RXM200 Massager is not a children's toy.

The RXM200 Massager is designed exclusively for massaging body areas. Store your RXM200 Massager in a clean and dry place protected from major temperature changes and other potentially harmful factors. Additionally, the RXM200 Massager is designed for indoor use only.

- Do not use the RXM200 Massager when it is wet or damp.
- Dismantling or rebuilding the RXM200 Massager may cause damage and injury.
- Use the RXM200 Massager for approx. 15 minutes and then switch it off for approx. 20 minutes.
- If the power adapter or the adapter is damaged, the RXM200 Massager must not be used!

The manufacturer or dealer assumes no liability for damage caused by improper or incorrect use.

WARNING

The following people should not use the RXM200 Massager:

- People with osteoporosis
- People with a pacemaker
- People with heart disease
- Pregnant women and people with wounds and/or injuries
- And children under 12 years

FUNCTION

The RXM200 Massager is deemed suitable for the following people:

- People sitting / standing in the same position for a long time (for example sitting at a desk in the office for a long time)
- Nervous and / or tense people
- People who suffer from sleep disorders
- People with a stiff neck or stiff shoulder parts
- People with headaches or dizziness
- People with 'sticky' fascia

Before plugging into the main electricity

1. Remove the packaging material and all protective films.
2. Unwind the cable of the power plug completely to straighten it.
3. Insert the power connector into the Holstecker plug of the power supply connector.
4. Plug the power plug into a properly installed power outlet.

You can use the RXM200 Massager as shown for the whole body:

NECK

SHOULDER

ABDOMEN

BACK

LEGS

ARMS

FEET

SERVICE

The control Panel:

How to operate the RXM200 Massager:

- Place the massager on the desired body zone.
- Put your hands in the arm loops and keep them at a comfortable point of tension for as long as you are using the RXM200 Massager on the neck, shoulder, back or hips. Alternatively, the fixable straps can be used.
- Press the On / Off button on the control panel and hold for 3 seconds. The RXM200 starts with the lowest massage option and the highest heat support (also when in battery mode).
- Use the button to increase the speed of the massage function. The blue glow of the button symbolizes the higher level of the massage and pump function (also when in battery mode).
- During the massage, you can adjust the heat with the button. The heating function is preset by default at the highest level. The white light of the button symbolizes the switched on heat function, in the lowest level. The blue glow of the key symbolizes the next higher level. The heat support is switched off as long as the button does not light up anymore (also when in battery mode).
- Use the button to turn ‚pump-mode‘ on and off. The pump mode is shown by default. The blue glow of the button symbolizes the deactivated pump assistance (also when in battery mode).

! TROUBLE-SHOOTING

- If the RXM200 Massager can not be switched on, please check whether the plug is correctly inserted and check the power supply.
- If the sound of the RXM200 Massager gets progressively louder, it may be because it is in continuous use for too long. Switch off the RXM200 Massager for at least 30 minutes before restarting.
- If the RXM200 Massager stops during use, it may be because the engine has overheated. Leave the RXM200 Massager turned off for at least 30 minutes before restarting.

Problem	Possible cause	Problem Resolution
The RXM200 Massager makes strange or loud noises.	The massage balls are electrically operated such that their activation is audible.	No correction necessary.
	The noise becomes increasingly louder.	Switch off the RXM200 Massager for at least 30 minutes.
	The massage balls rub on the fabric during normal operation.	No correction necessary.
The RXM200 Massager does not work.	The connectors are not connected properly.	Check the connectivity of the power connector and the power supply connector.
	The battery is not fully charged. Battery empty.	Fully charge the battery.

IMPORTANT: *In the event of a fault, do not attempt any repairs yourself. There are no warranty claims in such cases.*

WARRANTY

You have a 2-year warranty against defects in both material and workmanship, provided that the RXM200 Massager has been used as prescribed.

The warranty does not apply to damage caused by:

- Accident or unforeseen events (such as lightning, water, fire, etc.).
- Improper use or transportation.
- Disregard of the safety and maintenance instructions.
- Other improper uses, processing or modifications.

General information

The RXM200 Massager is intended for private use only and not for medical or commercial use. It is suitable for full body massage: neck, shoulder, back, legs, stomach, feet and arms. If you have any health concerns, please talk to your doctor before using the RXM200 Massager.

The massage should always be a pleasant experience. If you find the massage uncomfortable or causing pain, please change the intensity of the massage and/or adjust the tension of the arm loops to achieve a satisfactory level of comfort.

Do not use the RXM200 Massager while driving.

After 15 minutes of uninterrupted massage, the device switches off automatically.

Cleaning and Care

Before cleaning the RXM200 Massager, be sure to fully turn off the power and unplug the adapter.

Please note:

- Note that the RXM200 Massager is not washable per se!
- Always use a dry cloth for cleaning.
- Under no circumstances should any liquid get into the device or the accessories, otherwise a short circuit may occur.
- Minor, small stains may be removed only with a slightly damp cloth or sponge (with possibly a little soapy water).
- Do not use aggressive or solvent based cleaners. Do not use brushes with metal or nylon bristles, as well as sharp or metallic cleaning items such as knives and similar sharp objects.
- After cleaning, do not use the RXM200 Massager until it has completely dried.
- Keep the RXM200 Massager in a protective cover when not in use to prevent the accumulation of dust.

Technical data of the AC adapter

Manufacturer's name or trade mark, trade register number and address	ZHIMAI ELECTRONIC TECHNOLOGYLTD No.2 Fengyuan Road, Zhenglou Country, Wanquan Town, Wenzhou City, Zhejiang Province, China		
Model Identifier	CP1225		
	Measured value	Measuring tolerance	Unit
Input Voltage	230	+/- 1%	V
Input AC Frequency	60	+/- 1%	Hz
Output Voltage	11.93	+/- 1%	V
Output Current	2000	/	A
Output Power	23.86	/	W
Average active efficiency	86.08	+/- 2%	%
Efficiency at light load (10%)	83.65	/	%
No-load power consumption	/	/	W

Recycling loop

Packaging material can be conveyed back to the raw material cycle. Dispose of packaging material in accordance with legal provisions. Information can be retrieved from the return or collections systems of your community.

Notes on the Electrical and Electronic Equipment Act (ElektroG)

We would like to point out to owners of electrical and electronic equipment that, according to the applicable legal regulations, electrical devices must be disposed of separately from municipal waste. The symbol of a crossed-out waste bin shown on the electrical and electronic equipment shown also indicates the obligation to dispose of the batteries separately.

According to the ElektroG, as distributors of electrical devices, we are obliged to take back certain old electrical and electronic equipment free of charge in accordance with Section 17 (2) ElektroG and in conjunction with Section 17 (1) ElektroG. We comply with this legal obligation through DR Deutsche Recycling Service GmbH / take-e-way GmbH. You can also hand in your electrical equipment at the designated collection points of the public waste disposal authorities. By disposing of your old devices according to the legal requirements, you make an important contribution in the reuse and recycling of old devices.

As the end user, you are responsible for the deletion of all manner of personal data on electronic devices before handing the devices over for disposal.

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We / Wir, the / die Sportstech Brands Holding GmbH, 10785 Berlin, Germany
Address / Anschrift: Sportstech Brands Holding GmbH, Potsdamer Platz 11, D-10785, Berlin, Germany
Sporttech Brands Holding GmbH